

## Dominant Flex Surgical Suction Pump

### POWERFUL AND RELIABLE SUCTION



#### Powerful

- High-power flow rate of up to 60l/min while operating whisper-quiet at <math><40\text{dB(A)}^\*</math>
- Constant vacuum – no fluctuations
- Superfast ramp-up

#### Reliable

- Extreme durability
- Swiss made, high-quality components

**Independent, powerful and reliable – a strong alternative to wall vacuum**

The Dominant Flex serves nearly every suctioning need in the hospital, clinic or doctor’s practice. It fulfills the requirements for cardiac surgery thanks to Cardiac Floating rating. The pump offers a strong alternative to wall vacuum, as the independent vacuum avoids overflow into wall systems and provides a smooth, consistent, and reliable vacuum level. The heavy-duty piston-cylinder system creates whisper-quiet and robust operation. The initial investment in this suction pump is easily recovered from the many years of reliant functioning with minimal maintenance requirements and long replacement cycles.

The Dominant Flex offers flexibility, allowing the user to select a high, mid or low flow rate and thereby adjust the vacuum build-up time and sound level. High flow and high vacuum power also make it the perfect liposuction aspirator.

**Key features at a glance**



**Powerful but silent**

- ④ Powerful drive unit is whisper quiet at < 40 dB(A)\*.
- ⑤ Patented membrane regulator creates fast vacuum build-up.
- ⑥ Dominant Flex gives you the freedom to choose the required flow.





**User-centric design for ease of use**

- ⑦ Mechanical overflow set protects the pump.
- ⑧ Clean Touch elements and the housing with no gaps or grooves facilitate easy cleaning and minimize risk of infections.
- ⑨ Integrated foot on/off switch enables hands free operation.



**Designed to last a lifetime**

- ① Heavy-duty duty piston/cylinder technology offers exceptional durability.
- ② With 50 years of experience, Medela has created robust pumps with a 5-year warranty on the complete pump.
- ③ Patented flat springs minimize mechanical wear.

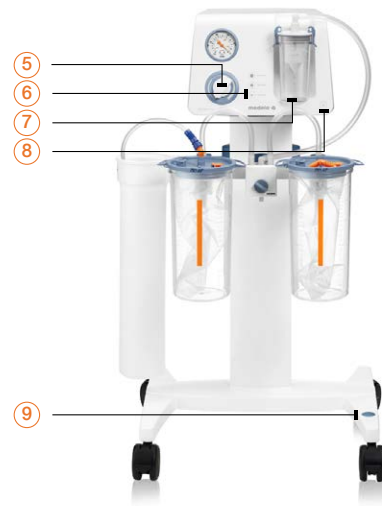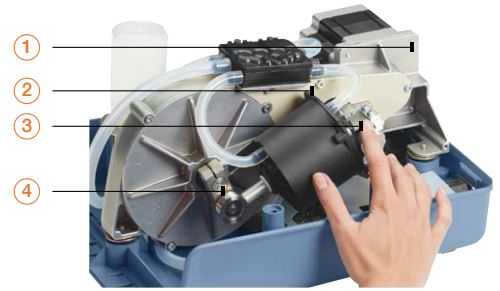

For an entry-level solution on surgical suction, **Medela’s Basic suction pump** offers ease of use and durability thanks to simplified, unique design (High vacuum – 90 kPa/ 675 mmHg; High flow: 30l/min.).

For increased mobility, **trolley** or portable versions are available. Medela’s fluid collection systems complement the pump enabling safe and silent collection of bodily fluids.

\*Data on file at Medela AG

**Technical data Dominant Flex**



high vacuum  
- 95 kPa/ 713 mmHg



9,3 kg (rack version)  
20.5 lbs

5 years warranty



high flow  
40/50/60 l/min.



H x W x D  
(rack version)  
210 x 305 x 375 mm

Medical Vacuum Technology  
for Healthcare Professionals

**Please contact us or your local  
Medela Representative for details.**

 Medela AG  
Lättichstrasse 4b  
6341 Baar, Switzerland  
www.medela.com

Local contact

